

INDIVIDUEL PRO

by Individuel Genève

USER GUIDELINE

Individuel
PRO
GENÈVE


INDEX

Step 1 Unboxing

Step 2 Assemble Technology Head & Capsule

Step 3 Protocol: Capsule treatment

Step 4: After Usage Instructions

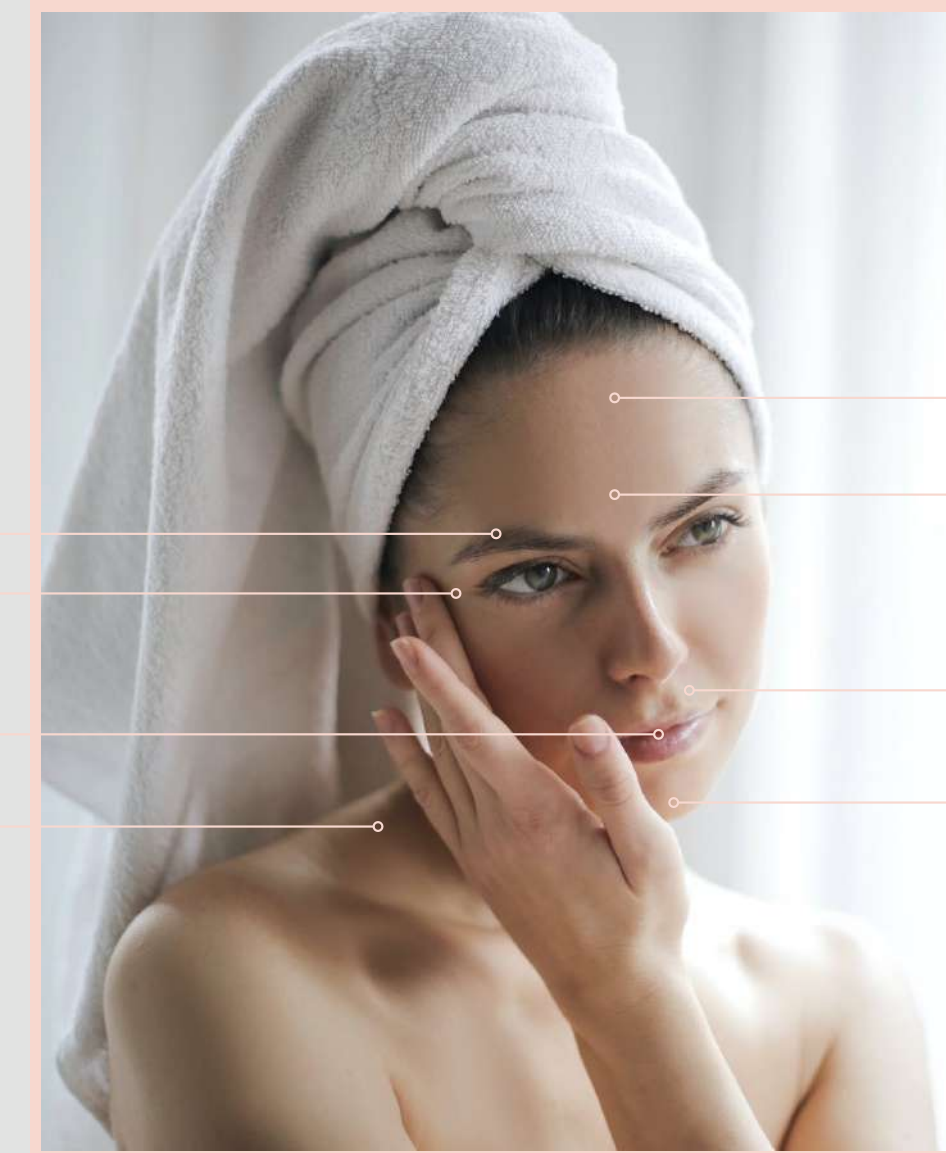

*Brow lift & asymmetric
eyebrow shaping*

*Crow's feet
(Peri-orbital lines)*

Thinning lips

*Neck folds
)Platysma Muscle(*

Forehead "worry" lines

Glabellar "Frown" lines

*Wrinkles around lips &
oral commissures*

Puckered Chin



STEP 1

UNBOXING

UNBOXING

Inside:

- 1 Pink suitcase
- 2 Handle device
- 3 Charging base
- 4 Charger & Cord
- 5 2 Technology Heads
- 6 5 capsules
- 7 Folder with warranty and guideline



STEP 2

ASSEMBLE TECHNOLOGY HEAD & CAPSULE



IMPORTANT!

Be sure that your handpiece is fully charged prior to initial use.

Keep the handpiece at the charger whenever it is not in use.

1


Connect Charger



2


Attach the technology head to the handpiece .

Gently **twist** the technology head clockwise until it clicks.

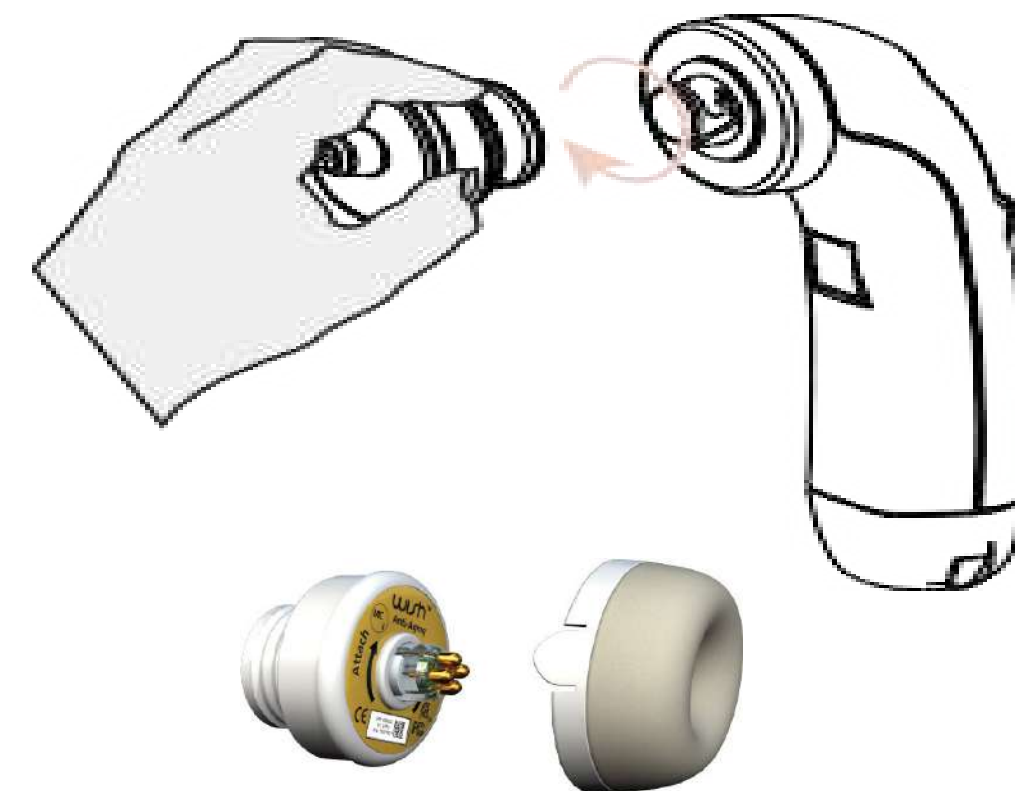

3

Hold the hand piece and **scan** the capsule against the technology head.

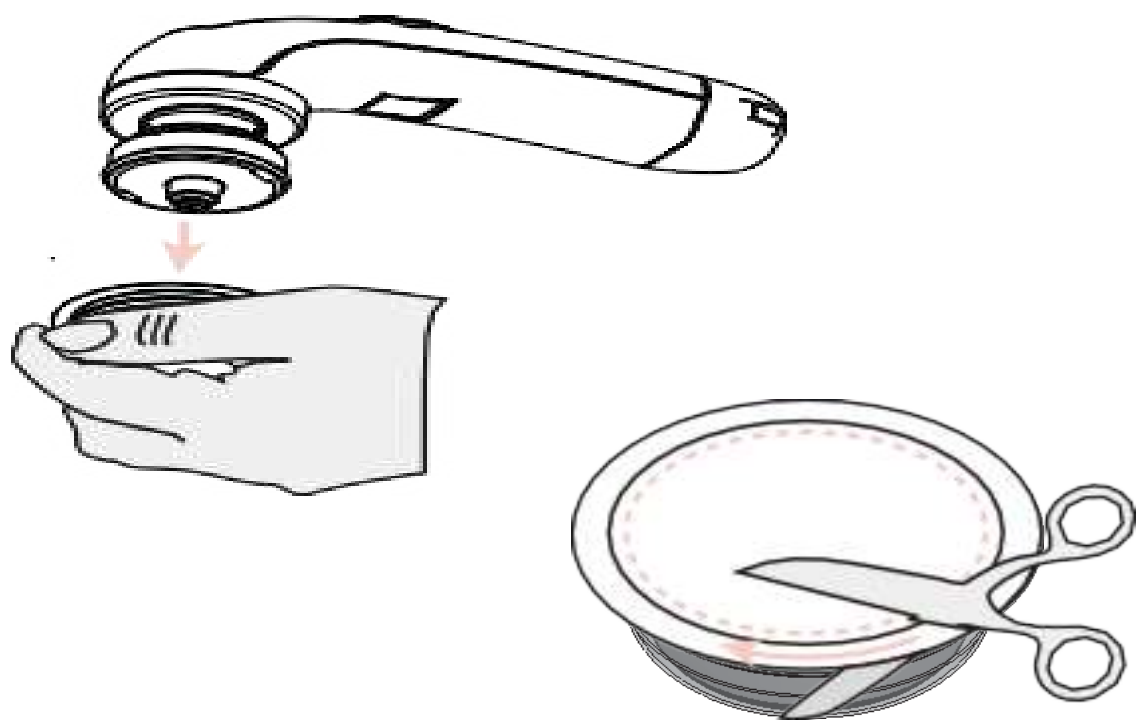

4

Open the appropriate capsule with a sharp tool .



5

Hold the capsules plastic and **attach** the capsule to the technology head.

Remove the plastic cup from capsule .
(Keep the plastic cup until the end).



If capsule shows VOID after scanning DO NOT open it. Contact us at: info@individuelcosmetics.com 0041 | 343 343 0800 | 78 86 810 22

5 Press the ON /OFF button to begin the procedure.



6 Adjust the intensity using the control buttons.
Press left arrow twice to reduce intensity to half circle.



PRESS POWER 
TWICE




PRESS LEFT ARROW 
TO REDUCE INTENSITY
TO HALF-CIRCLE



PRESS LEFT ARROW 
TO REDUCE INTENSITY
TO HALF-CIRCLE


STEP 3

PROTOCOL: CAPSULE TREATMENT

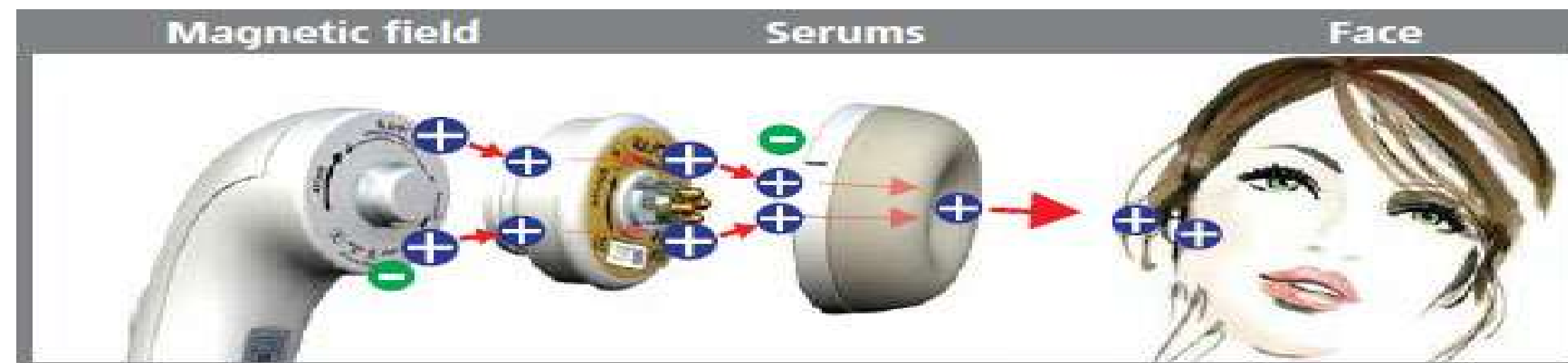

PROTOCOL: CAPSULE TREATMENT

- 1 Clean your face with face cleanser
- 2 Place the head and the selected capsule.
- 3 Slide over the face & neck & décolleté pressing gently, with an upwards and outwards motion for 8 minutes until full absorbency.
- 4 When finished, do not wash face for at least 2 hours after the treatment.

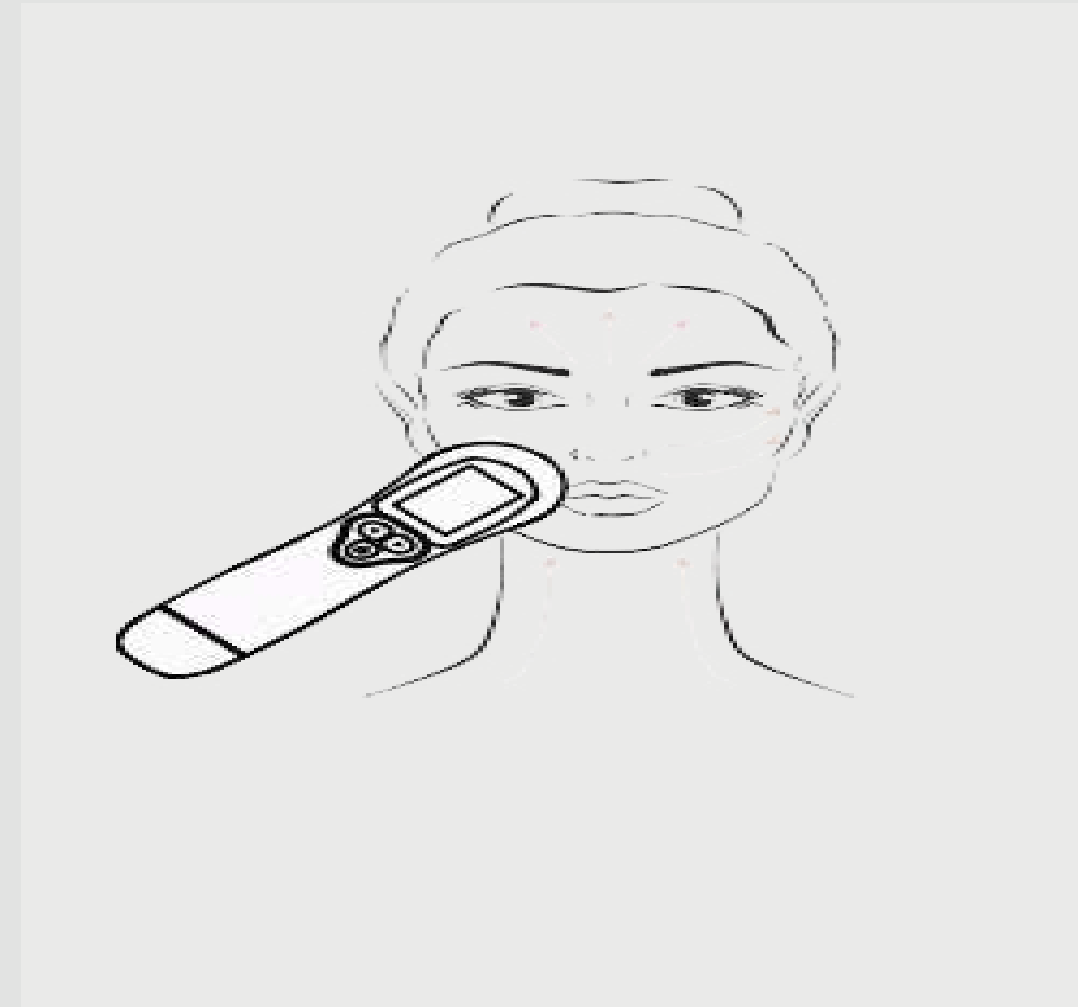



STEP 4


AFTER USAGE INSTRUCTIONS

IMPORTANT!


Do not grip the technology head by the rubber sleeve.

Be sure to grip the head on the hard plastic casing.

1 Shut off the device by pressing the ON/ OFF button.

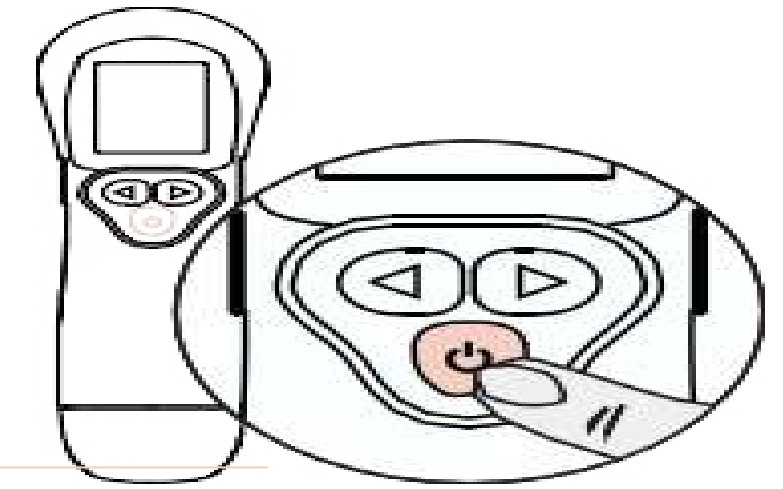


2 Firmly hold the handpiece in your right hand and remove the technology head with the capsule from hand piece by twisting it counter clock wise onto the plastic cup.

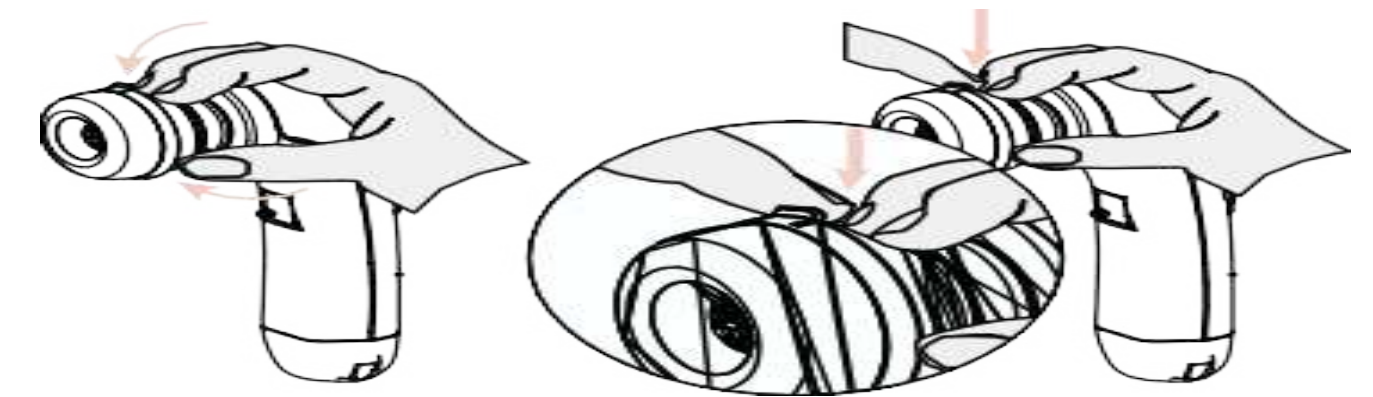


Insert the handpiece into the charger when not in use.

3

Detach the used capsule from the technology head and dispose it.

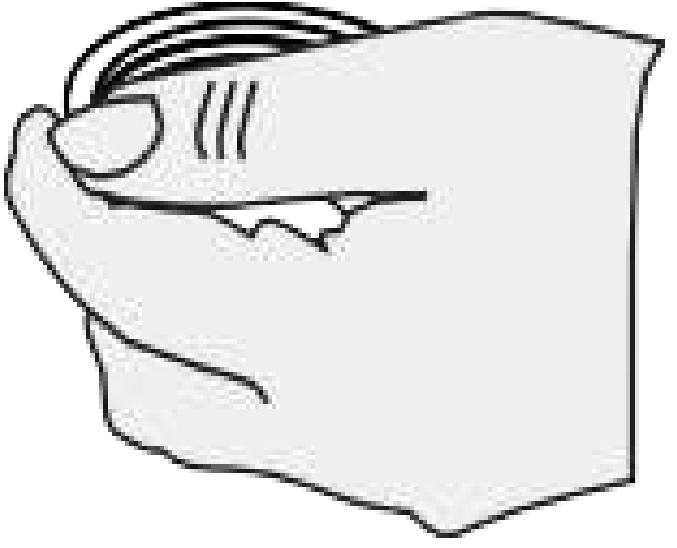

The capsules are designed for one- time use only.

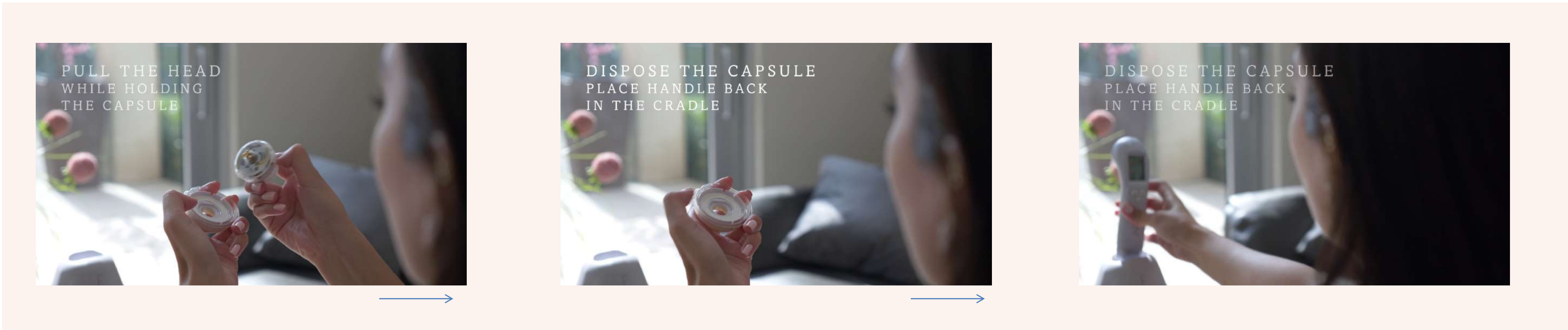

Insert the handpiece into the charger when not in use.

INDIVIDUEL PRO

Magnetic Field *Serum* *Face*
↓ ↓ ↓

